— B A C K S —

DEX BACKREST HARDWARE

QUICK START GUIDE

Tool Needed:

- 10mm wrench


- tape measure or level


Assemblies Included:

Two sets of fasteners provided. Reference Fastener Key below.

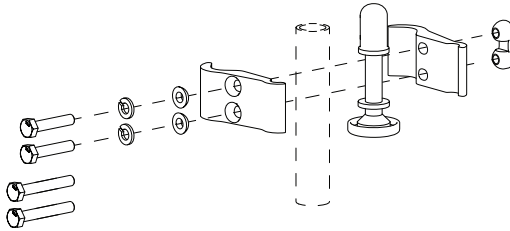
X2

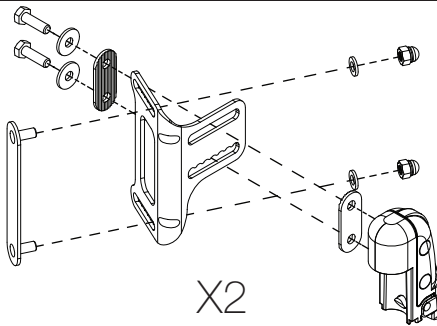
X2

Fastener Key:

Fastener drawings are actual size so hardware can be compared to drawings.

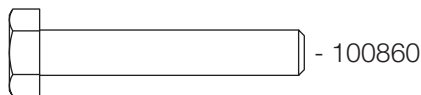
35mm (Used for 1" and 1-1/8" back canes)

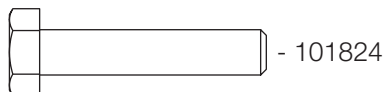30mm (Used for 3/4" and 7/8" back canes)


General Setup (Before Delivery)

- ① Determine installation slots based on configuration.

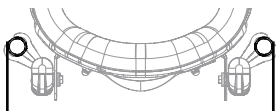
Fits Back Frame Widths		
Axiom Back Width	Inner Slots	Outer Slots
12" (30.5cm)	1" (28.0cm)	12" - 14" (30.5cm - 35.5cm)
14" (35.5cm)	12" - 13" (30.5cm - 33.0cm)	14" - 16" (35.5cm - 40.5cm)
16" (40.5cm)	14" - 15" (35.5cm - 38.0cm)	16" - 18" (40.5cm - 45.5cm)
18" (45.5cm)	16" - 17" (40.5cm - 43.0cm)	18" - 20" (45.5cm - 50.5cm)
20" (51.0cm)	18" - 19" (46.0cm - 48.5cm)	20" - 22" (50.5cm - 55.5cm)
22" (56.0cm)	20" - 21" (51.0cm - 53.5cm)	22" - 24" (56.0cm - 61.0cm)

- ② Remove backrest shell from cover. Install backrest bracket to shell with stud plate (A), two washers (B) and two dome nuts (C). Install so stud plates are 1" from the top of each bracket slot. Loosely tighten hardware. Repeat on opposite side ensuring both sides are set in the same Alignment Guide position.

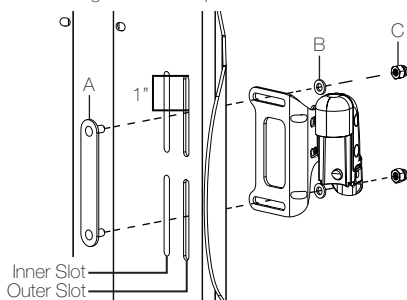
- ③ Hold the backrest shell against the chair at approximate height desired. Find a back cane location that can accommodate the clamp near the center of the slots on the backrest bracket. Typically the backrest will rest 1-2" above the cushion.

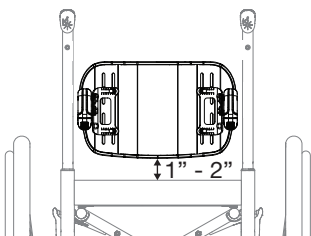
- ④ Determine the clamp fasteners needed by referencing the Fastener Key on the previous page of this document.


- ⑤ Loosely assemble the clamps onto the back canes at the desired location talked about in step 3. Bolts should thread inwards from the outside of the chair. Tighten clamp enough so it stays in place and repeat on opposite side.

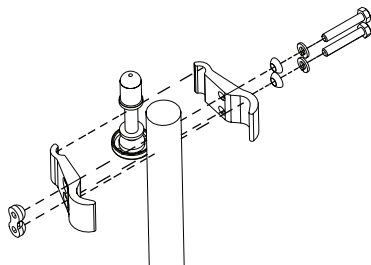
- ⑥ Use a level and a tape measure to ensure the clamps are level and the distance from backcane to clamp mount is 1/2".

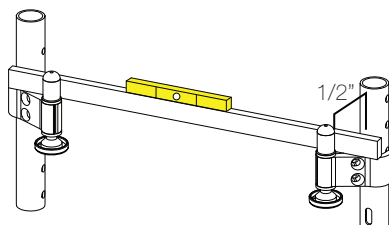- ⑦ Install backrest assembly onto the clamp mounts. Ensure both backrest plates are set in the same Alignment Guide position. Adjust the back height by moving the backrest along the plate up or down to desired height.

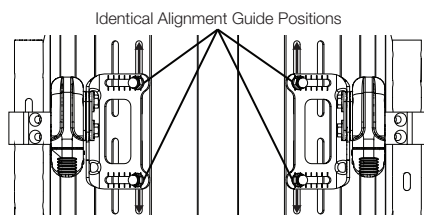- ⑧ Adjust the depth and angle by moving the backrest assembly to the desired positions. With position determined, ensure the bolts are secured in the same scallop on the bracket on both sides of the chair.

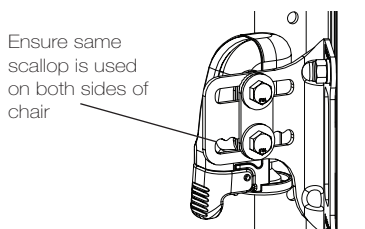
Fitting Adjustments (User Present)


⚠️ **WARNING** ⚠️

WARNING: IT MAY BE POSSIBLE TO ADJUST THE FINAL POSITION OF THE AXIOM BACK WITH AN END USER SEATED IN HIS WHEELCHAIR. THIS REQUIRES THAT THE END USER BE ABLE TO SIT IN A STABLE POSITION WITHOUT THE USE OF THE BACKREST. IT MAY ALSO BE POSSIBLE IF THERE ARE PEOPLE WHO CAN SUPPORT THE END USER DURING ADJUSTMENT. THIS SHOULD BE VERIFIED BY A CLINICIAN TRAINED TO MAKE THAT DETERMINATION. ALWAYS MAKE SURE THE WHEELCHAIR WHEEL LOCKS AND ANTI-TIPS ARE SET AND LOCKED. IF THESE CONDITIONS CANNOT BE MET, THE END USER SHOULD TRANSFER FROM THE WHEELCHAIR BEFORE MAKING ADJUSTMENTS. FAILURE TO TAKE THESE PRECAUTIONS MAY RESULT IN A FALL AND THE POTENTIAL FOR A SERIOUS INJURY.

- ① Adjust the angle and depth of the backrest by loosening the two bolts on the bracket. Grab and move/tilt the backrest assembly to desired position. Ensure both sides are set in the same scallop to ensure identical positioning on both sides. Retighten hardware once desired configuration is achieved.

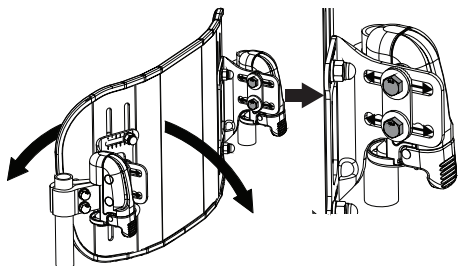- ② Adjust the height of the backrest by loosening the two bolts on the backrest bracket. Slide the shell up or down to desired height and ensure bolts are set in the same Alignment Guide position to ensure identical positioning on both sides. Retighten hardware once desired configuration is achieved.

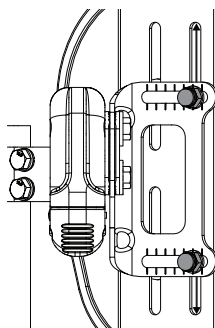Perform Check -

- Tighten all fasteners completely and assure that the hardware does not move out of the adjusted position under the load of the user.
- Verify the back releases and latches easily.
- Lift the backrest up by the strap and ensure it stays locked in.

If you have any questions or concerns about any aspect of this Axiom back, please do not hesitate to contact us by:

Telephone:

715-254-0991

In Writing at:

Ki Mobility

5201 Woodward Drive
Stevens Point, WI 54481
U.S.A

Or via email at:

sales@kimobility.com

Ki Mobility

5201 Woodward Drive
Stevens Point, Wisconsin 54481

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