— B A C K S —

BTX BACKREST HARDWARE

QUICK START GUIDE

Tool Needed:

- 10mm wrench


- 5mm Allen wrench


- tape measure or level


Clamp Positioning Chart

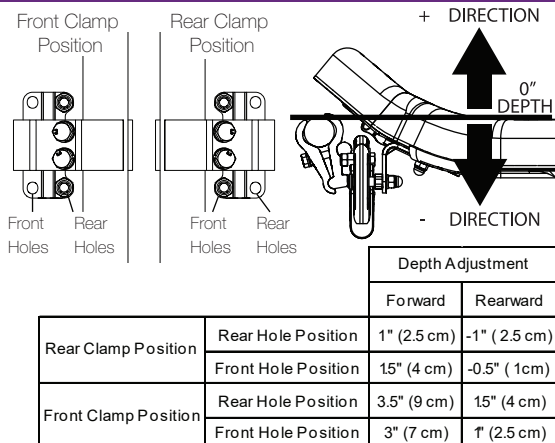Assemblies Included:

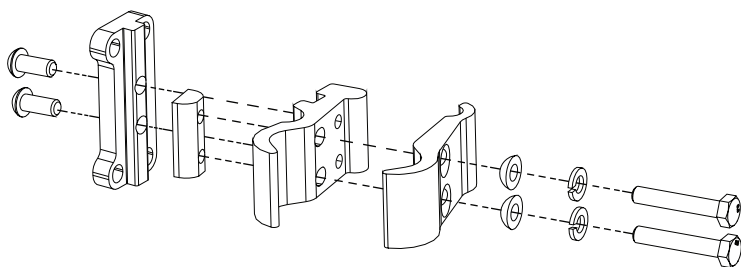
X2

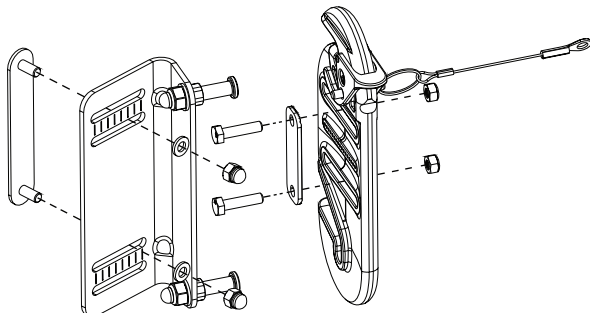X2


General Setup (Before Delivery)

- ① Determine installation slots based on configuration.

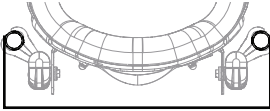
Fits Back Frame Widths		
Axiom Back Width	Inner Slots	Outer Slots
12" (30.5cm)	1" (28.0cm)	12" - 14" (30.5cm - 35.5cm)
14" (35.5cm)	12" - 13" (30.5cm - 33.0cm)	14" - 16" (35.5cm - 40.5cm)
16" (40.5cm)	14" - 15" (35.5cm - 38.0cm)	16" - 18" (40.5cm - 45.5cm)
18" (45.5cm)	16" - 17" (40.5cm - 43.0cm)	18" - 20" (45.5cm - 50.5cm)
20" (51.0cm)	18" - 19" (46.0cm - 48.5cm)	20" - 22" (50.5cm - 55.5cm)
22" (56.0cm)	20" - 21" (51.0cm - 53.5cm)	22" - 24" (56.0cm - 61.0cm)

- ② Remove backrest shell from cover. Install the backrest bracket to the shell with the stud plate (A), two washers (B) and two dome nuts (C). Install so stud plates are 1" from the top of each bracket slot. Loosely tighten hardware. Repeat on opposite side. Note: Ensure lanyard is at the top of the plate (Not seen below as it is mounted on the inside of the bracket on the top connecting rod hardware).

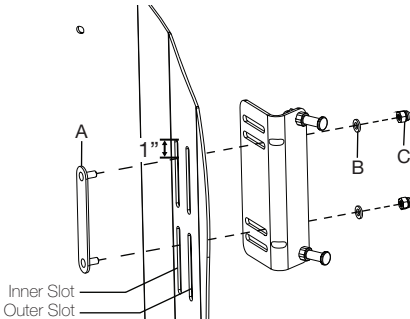
- ③ Hold the backrest shell against the chair at approximate height desired. Find a back cane location that can accommodate the clamp near the center of the slots on the backrest bracket. Typically the backrest will rest 1-2" above the cushion.

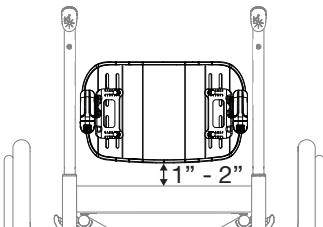
- ④ Determine the clamp orientation, forward or rear facing and the clamp bracket (Step 5 Image:A) orientation and hole use (See Clamp Positioning Chart on the second page of this document).

- ⑤ Loosely assemble the clamps onto the back canes at the desired location talked about in step 3. Bolts should thread inwards from the outside of the chair. Tighten clamp enough so it stays in place and repeat on opposite side.

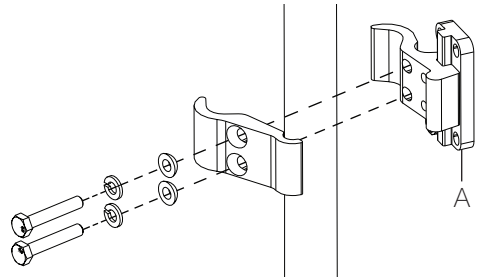- ⑥ Assume that the assemblies on each side are set at the same height and that the clamp brackets (Step 5 Image:A) are parallel to each other and set in the same orientation.

- ⑦ Install backrest mount to clamp assembly with two bolts (C), saddle (B) and two nuts (A) using a 10mm wrench. Repeat on opposite side.

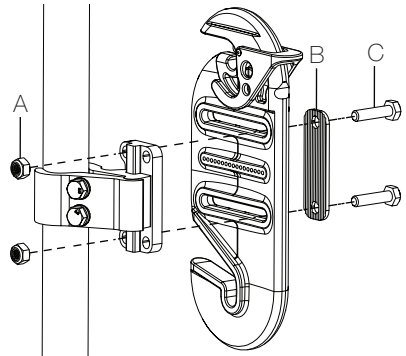
- ⑧ Install the backrest to the clamp and backrest mount assembly.

- ⑨ Adjust the backrest height to desired height by moving the backrest up or down along the backrest bracket hardware. See Step 2 of custom adjustments for imagery.

- ⑩ Adjust the backrest depth and angle to desired position by moving the backrest uniformly along the hardware attaching the clamp to the receiver on both sides. See Step 1 of Custom Adjustments for imagery.

- ⑪ Install the lock pin.

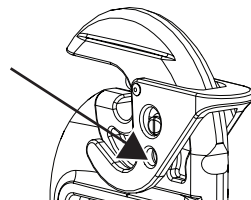Perform Check, see Check Section in next section.

Fitting Adjustments (User Present)



WARNING

WARNING: IT MAY BE POSSIBLE TO ADJUST THE FINAL POSITION OF THE AXIOM BACK WITH AN END USER SEATED IN HIS WHEELCHAIR. THIS REQUIRES THAT THE END USER BE ABLE TO SIT IN A STABLE POSITION WITHOUT THE USE OF THE BACKREST. IT MAY ALSO BE POSSIBLE IF THERE ARE PEOPLE WHO CAN SUPPORT THE END USER DURING ADJUSTMENT. THIS SHOULD BE VERIFIED BY A CLINICIAN TRAINED TO MAKE THAT DETERMINATION. ALWAYS MAKE SURE THE WHEELCHAIR WHEEL LOCKS AND ANTI-TIPS ARE SET AND LOCKED. IF THESE CONDITIONS CANNOT BE MET, THE END USER SHOULD TRANSFER FROM THE WHEELCHAIR BEFORE MAKING ADJUSTMENTS. FAILURE TO TAKE THESE PRECAUTIONS MAY RESULT IN A FALL AND THE POTENTIAL FOR A SERIOUS INJURY.

- ① To adjust the angle of the backrest, loosen the hardware attaching the clamp to the receiver on both sides of backrest. Tilt the backrest shell forward or backward to ensure uniform tilting on both sides. Tighten hardware once desired configuration is achieved.

- ② To adjust the backrest up or down, loosen the backrest bracket hardware (A) and move the assembly up or down. If more range is needed, the hardware can be removed and reinstalled into the other backrest bracket slots. Ensure both sides are set in the identical alignment guide positions in the backrest bracket. Retighten hardware once desired configuration is achieved.

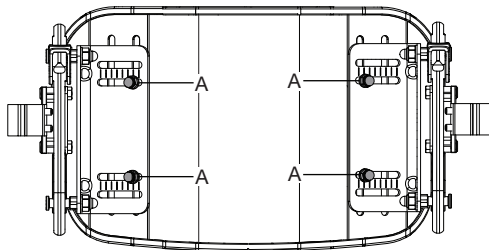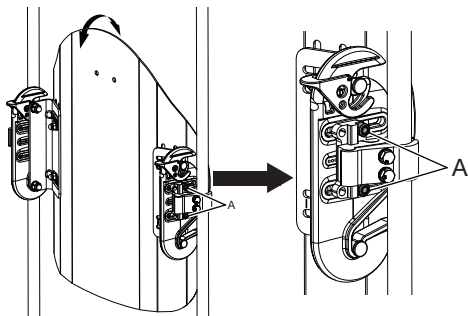Perform Check -

- Tighten all fasteners completely and assure that the hardware does not move out of the adjusted position under the load of the user.
- Verify the back releases and latches easily.
- Lift the backrest up by the strap and ensure it stays locked in.

If you have any questions or concerns about any aspect of this Axiom back, please do not hesitate to contact us by:

Telephone:

715-254-0991

In Writing at:

Ki Mobility

5201 Woodward Drive
Stevens Point, WI 54481
U.S.A

Or via email at:

sales@kimobility.com

Ki Mobility

5201 Woodward Drive
Stevens Point, Wisconsin 54481

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