— B A C K S —

FXT BACKREST HARDWARE

QUICK START GUIDE

Tool Needed:

- 10mm wrench

- 5mm Allen wrench


- tape measure or level


Assemblies Included:

Two sets of fasteners provided. Reference Fastener Key below.

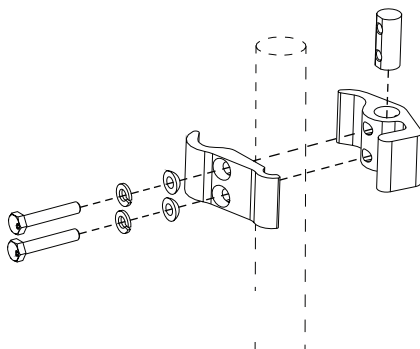
X2

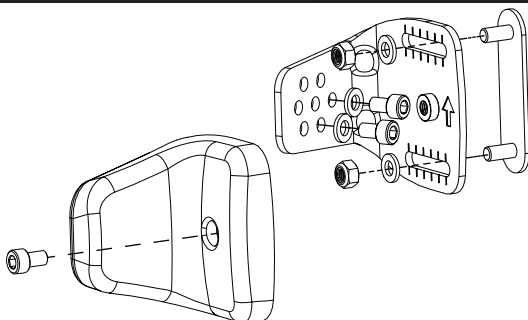X2


General Setup (Before Delivery)

- 1 Determine installation slots based on configuration.

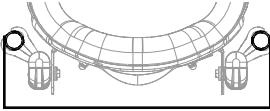
Fits Back Frame Widths		
Axiom Back Width	Inner Slots	Outer Slots
12" (30.5cm)	1" (28.0cm)	12" - 14" (30.5cm - 35.5cm)
14" (35.5cm)	12" - 13" (30.5cm - 33.0cm)	14" - 16" (35.5cm - 40.5cm)
16" (40.5cm)	14" - 15" (35.5cm - 38.0cm)	16" - 18" (40.5cm - 45.5cm)
18" (45.5cm)	16" - 17" (40.5cm - 43.0cm)	18" - 20" (45.5cm - 50.5cm)
20" (51.0cm)	18" - 19" (46.0cm - 48.5cm)	20" - 22" (50.5cm - 55.5cm)
22" (56.0cm)	20" - 21" (51.0cm - 53.5cm)	22" - 24" (56.0cm - 61.0cm)


- 2 Remove backrest shell from cover. Install backrest bracket to shell with stud plate (A), two washers (B) and two nuts (C). Stud plates should be 1.5 - 2" from the top of the slot. Loosely tighten hardware. Repeat on opposite side. Note: For FXT Hardware only, the stud plate can be installed so both studs are in the same slot if the configuration requires.

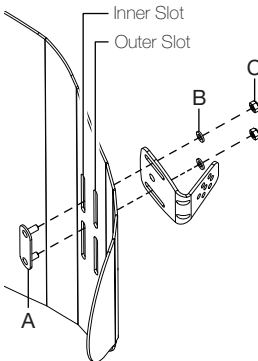- 3 Hold the backrest shell against the chair at approximate height desired. Find a back cane location that can accommodate the clamp near the center of the slots on the backrest bracket. Typically the backrest will rest 1-2" above the cushion.

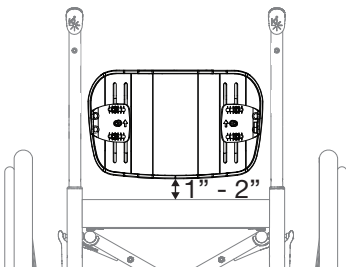
- 4 Choose the desired angle and depth set up. Use the images and charts below to help determine how to achieve desired configurations.

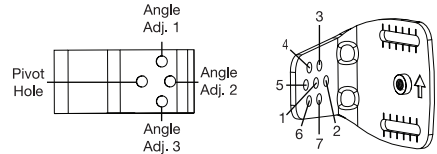	-0°		-2.5°		-5°	
Bracket Arrow ↑	Pivot	Angle 3	Pivot	Angle 2	Pivot	Angle 1
Rear Position (Closer to back canes)	1	6	1	5	1	4
Front Position (Farther forward from back canes)	2	7	2	1	2	3
	+0°		+2.5°		+5°	
Bracket Arrow ↓	Pivot	Angle 3	Pivot	Angle 2	Pivot	Angle 1
Rear Position (Closer to back canes)	1	6	1	5	1	4
Front Position (Farther forward from back canes)	2	7	2	1	2	3

- 5 Install backrest clamp assembly to the backrest bracket with two screws (A) and two washers (B) using a 5mm Allen wrench. Note: The angle is set by installing the hardware in the holes determined in Step 4.

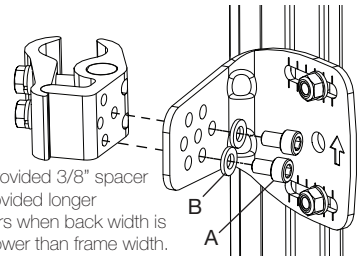
*Use provided 3/8" spacer and provided longer fasteners when back width is 2" narrower than frame width.

- 6 Adjust the depth and angle by moving the backrest assembly to the desired positions. With position determined, ensure the bolts are secured in the same scallop on the bracket on both sides of the chair.

✓ Perform Check, see Check Section in next section.

- 7 Loosely install the back mount cover with one screw using a 5mm Allen wrench and repeat on opposite side.

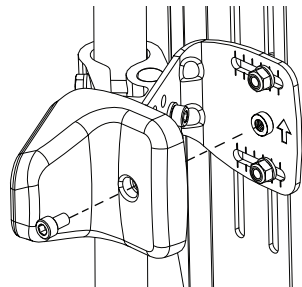
Fitting Adjustments (User Present)


⚠️ WARNING ⚠️

WARNING: IT MAY BE POSSIBLE TO ADJUST THE FINAL POSITION OF THE AXIOM BACK WITH AN END USER SEATED IN HIS WHEELCHAIR. THIS REQUIRES THAT THE END USER BE ABLE TO SIT IN A STABLE POSITION WITHOUT THE USE OF THE BACKREST. IT MAY ALSO BE POSSIBLE IF THERE ARE PEOPLE WHO CAN SUPPORT THE END USER DURING ADJUSTMENT. THIS SHOULD BE VERIFIED BY A CLINICIAN TRAINED TO MAKE THAT DETERMINATION. ALWAYS MAKE SURE THE WHEELCHAIR WHEEL LOCKS AND ANTI-TIPS ARE SET AND LOCKED. IF THESE CONDITIONS CANNOT BE MET, THE END USER SHOULD TRANSFER FROM THE WHEELCHAIR BEFORE MAKING ADJUSTMENTS. FAILURE TO TAKE THESE PRECAUTIONS MAY RESULT IN A FALL AND THE POTENTIAL FOR A SERIOUS INJURY.

- ① To adjust the angle and depth of the backrest, loosen the pivot screw on the backrest bracket and remove the angle adjustment screw. Grab and move/tilt the backrest assembly to desired position. Align holes for desired backrest angle and reinstall the angle adjustment screw. Retighten the pivot screw. Use the images and charts in General Setup Step 4 to identify screws and figure out the angle configurations.

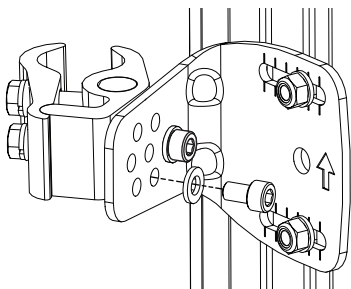- ② Remove the back mount cover. Loosen the backrest bracket hardware (A) to move the entire back assembly up or down in the slots to desired height. If more adjustability is needed, the hardware can be removed and reinstalled so both studs are in the same stud slot. Ensure both sides are set in the identical height configuration.

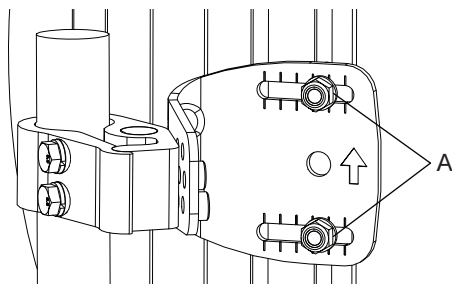- ✓ Perform Check -
a. Tighten all fasteners completely and assure that the hardware does not move out of the adjusted position under the load of the user.

- ③ Install the back mount cover with one screw using a 5mm Allen wrench and repeat on opposite side. See General Setup Step 7 for imagery if needed.

If you have any questions or concerns about any aspect of this Axiom back, please do not hesitate to contact us by:

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U.S.A