

SPAN PressureGuard® Bariatric APM



Bariatric Alternating Pressure Mattress System

Alternating Pressure therapy with gentle, effective pressure changes for user up to 750 lbs

Uniquely equipped to meet the extreme skin integrity challenges of bariatric care

Safety Edge™ for added fall prevention and stability

science, safety, simplicity


 **savaria**®

SPAN PressureGuard® Bariatric APM


Segmented top reduces sharp pressure drop-offs, while air system below handles the “heavy lifting,” alternating weight-bearing across the surface in user-selectable intervals. “Two up, one down” inflation pattern ensures maximum patient support throughout the cycle.

Treatment Flexibility

Two selectable therapy modes:

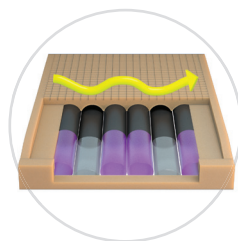


Alternating Pressure
Gently changes loading across the surface (10-minute cycle).

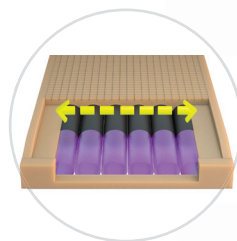

“Float” Mode
Provides active, powered immersion therapy.

PressureGuard Design

Integrated air system/engineered foam shell with “Safety Edge” bolster system delivers superior edge-of-bed safety, stability and comfort.

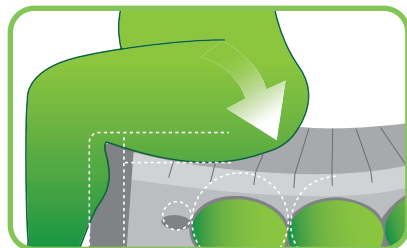

Safety Edge

Ensures safer transfers while minimizing the risk of falls and patient entrapment.



Unique Geo-Matt® Design

Creates tissue-friendly, anti-shearing surface with more than 800 individually responsive cells.



Digital Control Unit

A robust, sophisticated control unit provides a high level of customization to meet each user's unique needs. Features include:

- “Float” (powered immersion) mode
- Timed (20-minute) Auto Firm mode
- Ten electronically-controlled comfort settings
- Programmable cycle times of 20, 25 or 30 minutes
- Low pressure indicator, power failure indicator and audible alarm
- Lockout mode



Indications

- Aggressive treatment (through Stage 4)
- Prevention of pressure injuries in immobile or high-risk patients
- Can't or won't reposition frequently
- Can be positioned off the pressure injury



Durable Fluid-Resistant Cover

The bacteriostatic top fabric is flame-resistant, fluid-resistant, tear-resistant, and replaceable. Wipes clean in place with standard hospital-grade cleaners. The pleated design allows full integration with the mattress's shear-relieving surface while minimizing hammocking.

ITEM #	DESCRIPTION
AP8042-29 Bariatric APM	80"L x 42"W x 7"H
AP8048-29 Bariatric APM	80"L x 48"W x 7"H
AP8054-29 Bariatric APM	80"L x 54"W x 7"H
AP8442-29 Bariatric APM	84"L x 42"W x 7"H
AP8448-29 Bariatric APM	84"L x 48"W x 7"H
AP8454-29 Bariatric APM	84"L x 54"W x 7"H

CONTROL UNIT

Weight	6.3 lbs.
Voltage	110-120 AC
Frequency	60 Hertz
Dimensions	13.4" x 5.9" x 9.2"
Max Current	1.0 Amp
UL Listing	Certified to IEC 60601-1

WEIGHT LIMIT

120 lbs. - 750 lbs

MATTRESS WEIGHT

approx. 25 lbs.

WARRANTY*

Mattress	3 years
Cover	2 years
Control Unit	18 months

FLAMMABILITY

Conforms to NFPA 101 (Life Safety Code); Cal. TB 117; ASTM E1590; and 16 CFR Parts 1632 and 1633.



Get more information

800.888.6752

savariapatientcare.com



Span-America
Handicare USA

800.888.6752
866.891.6502

PO Box 5231 Greenville SC 29606
10888 Metro Court St. Louis MO 63043

Span-America is a division of Savaria Corporation, a global leader in accessibility. Learn more about Savaria at savaria.com.

PLspc-36 CC-5589 V. 1